


SAP Analysis Instructions



1


Open **HRx Plant Sap Analysis Sample kit** and verify following contents.

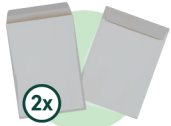

2x
Insulated envelopes



2x
Sap sample labels

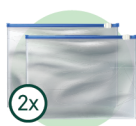

2x
Ziplock bags



2x
Ice packs



Sample Submission Form



Overnight label

2

You will also need:



Scissors



Scales

Take Sap Sample Pair:



NEW
Fully developed in the Top 30% of plant

Old
From viable bottom 30% of plant

*Be sure to include leaf and petiole

3


4



Take samples **1-3 hours from lights on** making sure that future samples of the same crop are taken at the same time.

6

Tare **bag and Weigh sample** to ensure 80 -100 grams of sample

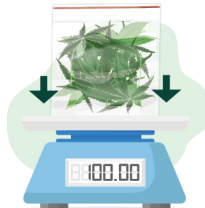

5




Place **NEW & OLD leaves** in separate zip-lock bags

*Be sure to let air out of bags after sampling

7

Place zip-loc inside **insulated packaging** and place properly filled out label onto appropriate sample

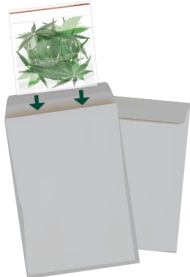

*Keep ice packs on outside of insulated bags

9

Package everything and ship same day as sampling overnight Mon-Wed.



8

Record appropriate **information** onto sample submission form

